



waterpik® waterflosser®

QUICK START GUIDE WP-150 Series


1 First Use

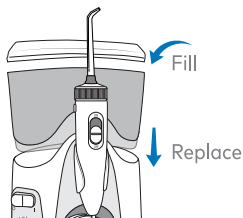

Fill reservoir with warm water and replace it firmly on the base.

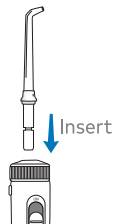

Insert tip firmly until it clicks into place.

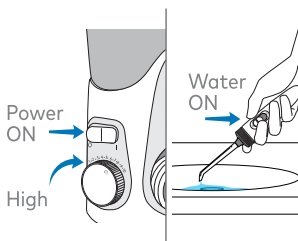


Turn pressure to high, point tip into sink. Turn Power switch ON and water switch ON until water flows. Turn OFF.


waterpik® waterflosser®

QUICK START GUIDE WP-150 Series


2 Water Floss

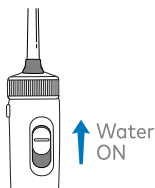

With Power OFF and water switch in OFF position, place tip in mouth. Start on low pressure and increase to your liking. Turn power ON. Turn water switch to ON.

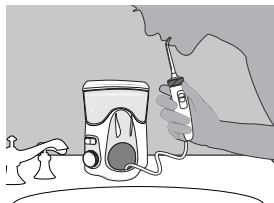

Lean over sink. Close lips enough to prevent splashing, while still allowing water to flow from mouth into sink.

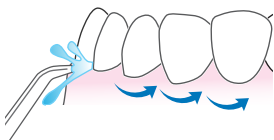

Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.