**waterpik®**


# **Waterpik® Sensonic® Professional Plus Toothbrush SR-3000**

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**Cepillo dental  
Waterpik® Sensonic®  
Professional Plus SR-3000**

**Spazzolino  
Waterpik® Sensonic®  
Professional Plus SR-3000**

**Зубная щетка Waterpik®  
Sensonic®  
Professional Plus SR-3000**

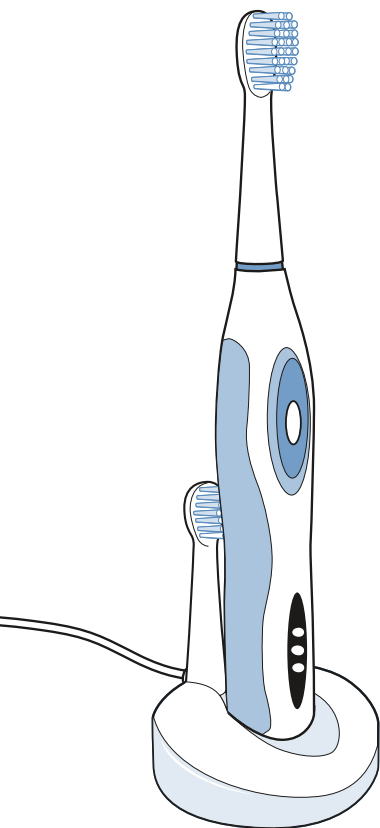.....  
[www.waterpik.com](http://www.waterpik.com)

## WELCOME

Congratulations, as an owner of the Waterpik® Sensonic® Professional Plus Toothbrush, you are now on your way to better plaque removal, naturally whiter teeth, and a healthier smile!

The Sensonic® Professional Plus Toothbrush utilizes an ultra high bristle speed and wide range of brushing motion. As a result, the Sensonic® Professional Plus Toothbrush generates superior cleaning action and delivers improved plaque removal, healthier gums, and a fresh clean feeling.

Thank you for choosing the Waterpik® Sensonic® Professional Plus Toothbrush.

Important Safeguards	4
Product Description	6
Getting Started	7
Care of your Waterpik® Sensonic® Professional Plus Toothbrush	10
Limited Two-Year Warranty	13

## IMPORTANT SAFEGUARDS

### IMPORTANT SAFEGUARDS


As with any electrical product, especially when children are present, always follow these basic safety precautions.

**PLEASE READ ALL INSTRUCTIONS BEFORE USING.**

### DANGER:

**To reduce the risk of electrocution:**

- Do not handle the plug with wet hands.
- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not immerse charger in water or other liquid.

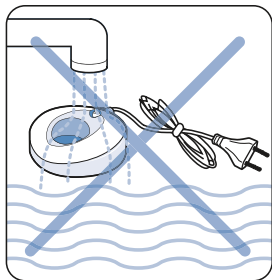- Do not handle charger base while bathing.
- Do not reach for a charger that has fallen into any liquid. Unplug immediately.
- Check the charger cord for damage before the first use and during the life of the product.

### WARNING:


**To reduce the risk of burns, electric shock, fire, serious injury or damage to tissue:**

- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Do not use this product if the charger has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if it has dropped into any liquid. Contact Waterpik International, Inc. or your local distributor which can be found at [www.waterpik.com/intl](http://www.waterpik.com/intl).
- Be sure only the bristles of the toothbrush touch your teeth and gums when using the device.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- For household use only. Do not use outdoors.
- Do not use while wearing any oral jewelry. Remove jewelry prior to use.
- Do not use if you have an open wound on your tongue.
- Use only those attachments, including toothbrushes and interdental brushes, which are recommended by Water Pik, Inc.

- Keep charging unit and cord away from heated surfaces.
- Do not operate or use where oxygen or aerosol sprays are being used.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult with your dentist or medical professional before using this instrument or any other oral hygiene aid.
- This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device containing magnets or electronics, contact your physician or the device manufacturer prior to using the Sonic Toothbrush.
- This product does not contain any user serviceable parts in the handle or charger. Do not attempt to open either component except for battery disposal. See page 12, "Battery Removal."
- Closely instruct and supervise children and individuals with special needs in the proper use of this product to ensure safe use.
- Children should be supervised to ensure that they do not play with the appliance.


**SAVE THESE INSTRUCTIONS**

## PRODUCT DESCRIPTION
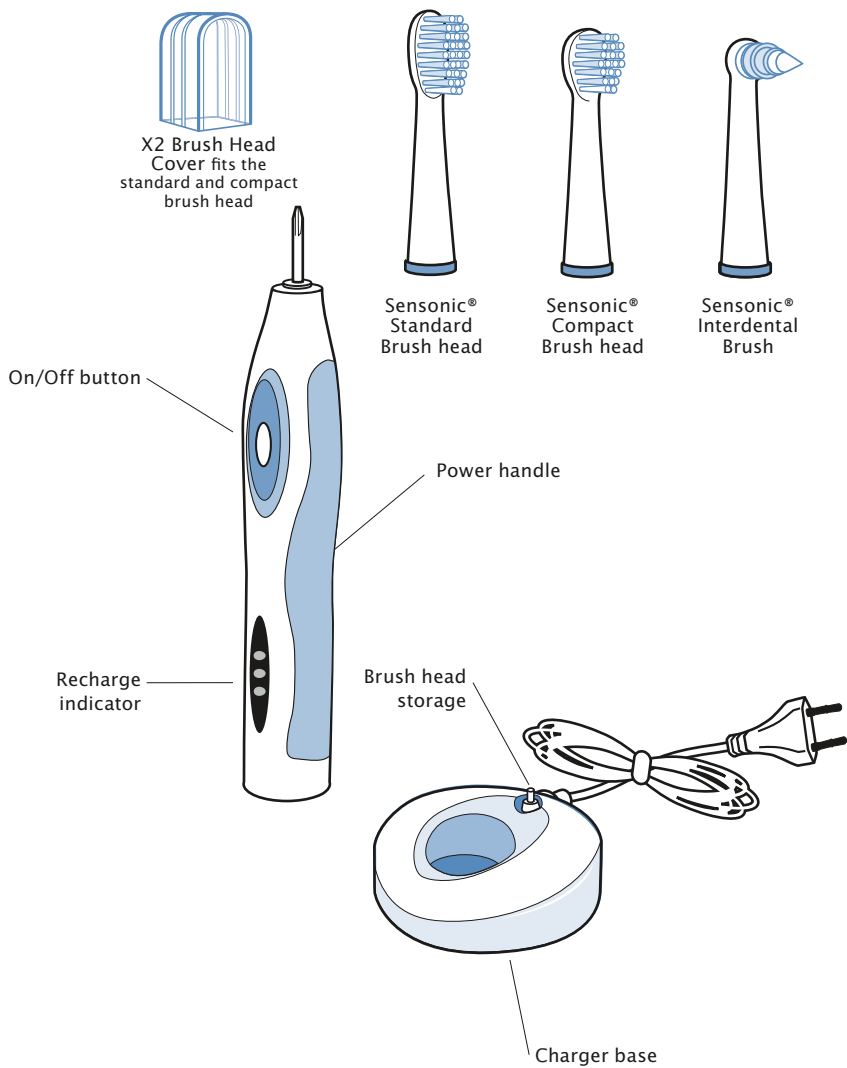
## Charging the Sonic® Professional Plus Toothbrush



**Charge for 24 hours before first use.** The charge indicator lights on the handle will go through a fill-up sequence when charging and all 3 lights will stay on continuously when the unit is fully charged. A full charge will last for approximately 1 week of typical use (2 minutes of brushing twice a day).

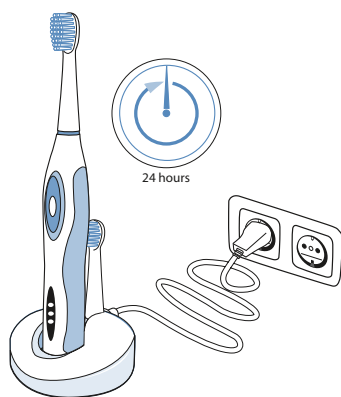If the toothbrush is to be stored for an extended period of time (six months or more), charge it for 24 hours before storing. Never store the toothbrush for extended periods unless it is fully charged.

## Recharge Indicator Lights

When the unit is taken off a powered charging base it will indicate the level of battery charge for 10 seconds and then the lights will go off. The lights will stay illuminated for 10 seconds after use indicating the level of battery charge.

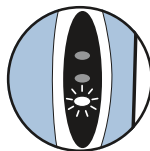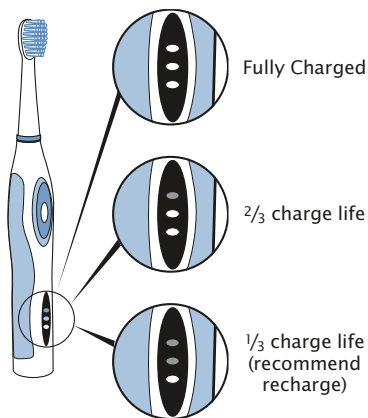Unit will not operate if the bottom LED begins to blink rapidly.

In order for the toothbrush to begin to work again, it must be placed on the charger for approximately 1 hour. It will still require a 24 hour charge to regain full charge status.


## GETTING STARTED

### Sonic brushing

- Slide the brush head onto handle and push it down firmly.


- Place the brush head into your mouth before turning the unit on.

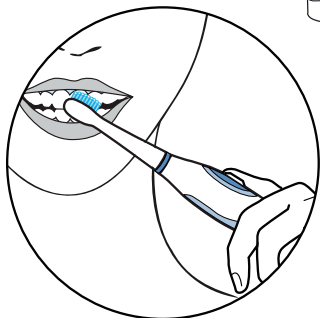- Select low or high brushing speed by pressing power button. Press button once for low speed, and twice for high speed.


- If you hold the on/off button down for approximately 1.5 seconds while in the low speed position the brush will shut off.
  - Low speed – gentle cleaning and gum massage
  - High speed – superior plaque removal
- For best results, position brush along gumline at a 45-degree angle.

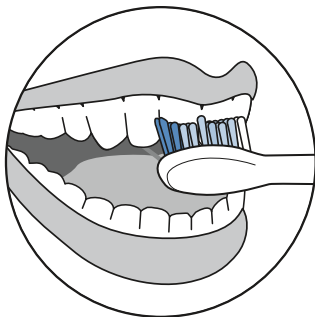
- Use a lighter brushing pressure than when using a manual toothbrush to optimize sonic cleaning effectiveness.

- Move the brush slowly around entire mouth. Many dentists recommend spending 30 seconds in each of the four quadrants of your mouth.
- Quad Timer: your Sensonic® Professional Plus Toothbrush will create a short pause every 30 seconds, to remind you to move to the next quadrant of your mouth.

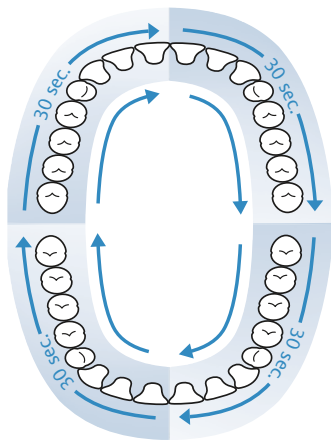
- 2-Minute timer: The brush will automatically shut off after 2 minutes of run time.

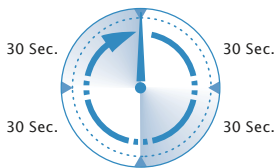
**2 minutes**

### Using the interdental brush

- Slide the brush head onto handle and push it down firmly.
- Place the brush head into your mouth before turning the unit on.
- Brush may be used with or without toothpaste.



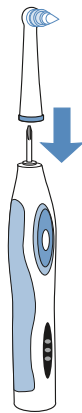
- Place brush tip at a 90-degree angle to the gumline between two teeth.
- Turn on power handle and lightly follow crevices between teeth.

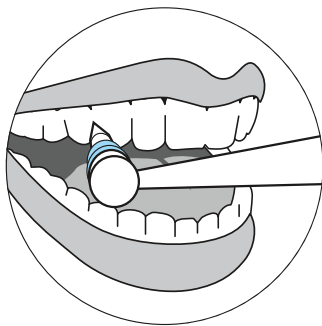
## Cleaning

1. After each use, turn off unit, remove the brush head from the handle and rinse brush head under running water.
2. Wipe the handle dry after each use.
3. Periodically, clean the charger base, first unplugging the unit and then wiping with a damp cloth.

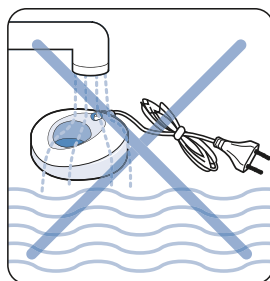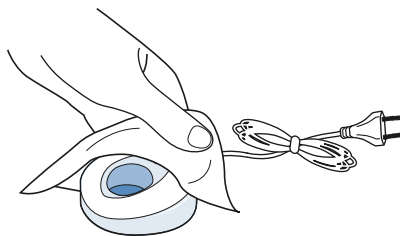
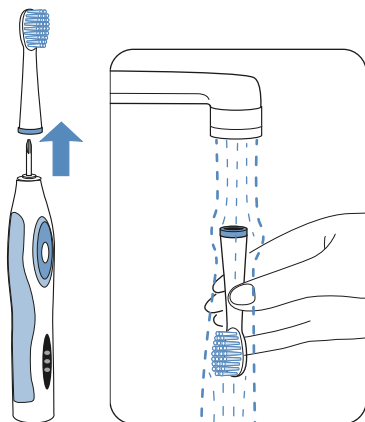
2. Wipe the handle dry after each use.
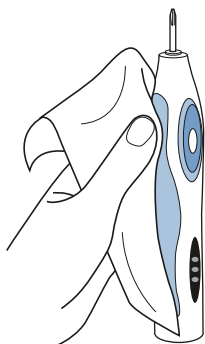
## Brush Head Cover

The brush head covers fit both the standard and small brush heads. For the standard brush head use the larger side and for the compact brush head use the smaller side.

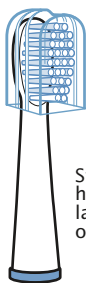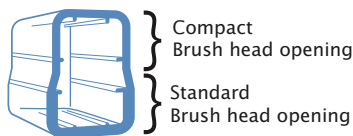Standard Brush head using the larger size opening

## Service

Waterpik® toothbrushes have no consumer serviceable or replaceable electrical items and do not require routine service maintenance.

For all your warranty and accessory needs, please visit [www.waterpik.com/intl](http://www.waterpik.com/intl) to find your local distributor. Please do not ship product to Waterpik International, Inc. offices, this will delay service. Refer to serial number and model number in all correspondence. These numbers are located on the bottom of your charger.

## Important Precautions

Using baking soda or peroxide containing toothpastes, abrasive or chemical cleaners may damage the plastic housing of this product.

## Disposal

This product contains rechargeable Nickel Metal Hydride batteries. The batteries are not replaceable.

Please do not dispose of the unit in the household waste. Remove the batteries per the instructions on page 12 and recycle at your local recycling center. In so doing, you will help preserve the environment.


To locate your local distributor, go to [www.waterpik.com/intl](http://www.waterpik.com/intl).

### Battery Removal **Note: This process is not reversible**


The rechargeable battery inside the toothbrush cannot be replaced but can be removed for recycling. Remove the battery only at the end of the product's life. Doing so earlier will destroy the product. Deplete the battery charge by running 2 minute cycles until the brush no longer operates. Observe basic safety precautions.

At the end of the product's useful life, open the handle as shown, remove the batteries and dispose of them according to local regulations. The rest of the product can be discarded per local regulations.


1. Use coin or washer to rotate cap counter clockwise until cap protrudes. Rotation will stop at this point.

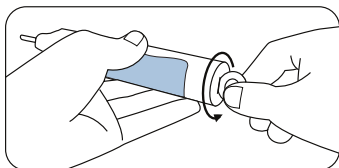
2. Push down hard with shaft on hard surface.

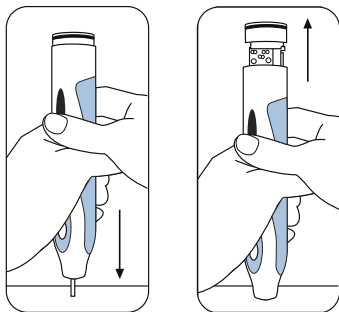3. Remove cap and pull internal body out. Cut tape and remove.

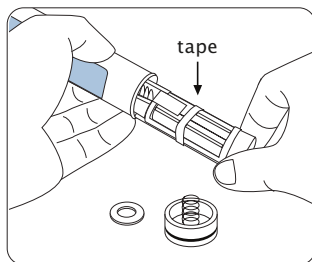4. Lift battery pack up and cut electrical wires, one at a time, to free from circuit board.
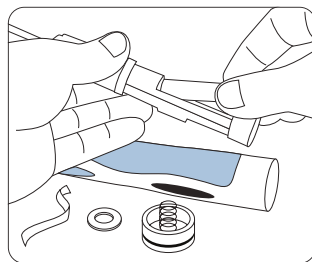
## LIMITED TWO-YEAR WARRANTY

Water Pik, Inc. warrants to the original purchaser/owner of this product that it is free from defects in materials and workmanship for two years from date of purchase. Save your receipt as your proof of purchase date. Water Pik, Inc. reserves the right to request the product to be returned for analysis to confirm the claim falls within our warranty policy. We will replace any part of the product, which in our opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase and was used according to these instructions. Repair or replacement services provided for by the guarantee will be carried out by the local distributor, whose details may be found on [www.waterpik.com/intl]. Upon request of the consumer, a copy of the guarantee complete with the local distributor's information will be provided by email or mail. This limited warranty excludes accessories or consumable parts such as tips, brush heads, etc.

This warranty gives you specific legal rights. You may also have other rights that vary depending on the law in your place of residence.

To locate your local distributor, go to [www.waterpik.com/intl](http://www.waterpik.com/intl).

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Waterpik® is a trademark of Water Pik, Inc. registered in Argentina, Australia, Benelux, Canada, Chile, China, Hong Kong, India, Japan, Korea, Mexico, New Zealand, Norway, Pakistan, Russian Federation, Singapore, South Africa, Switzerland, Taiwan, Ukraine, and the United States.

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Australia, Brazil, Canada, China, EU, Hong Kong, Japan, Mexico, Korea, Russian Federation, Switzerland, Trinidad & Tobago, and the United States.

Sonic® is a trademark of Water Pik, Inc. registered in Canada, China, EU, Japan, Switzerland, and the United States.

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©2012 Water Pik, Inc.
Sonic® è un marchio commerciale di Water Pik, Inc., depositato in Canada, Cina, UE, Giappone, Svizzera e USA.